


To avoid overheating of the stator coils the Optimus allows only 4 preparations within one minute. The 5th attempt will not lead to an acceleration. Ready disappears from the desk, but nothing happens. In addition 00XT+10LA will be logged to indicate e.g. a bouncing prep switch or "a nervous thumb".

Such symptoms can be avoided by programming

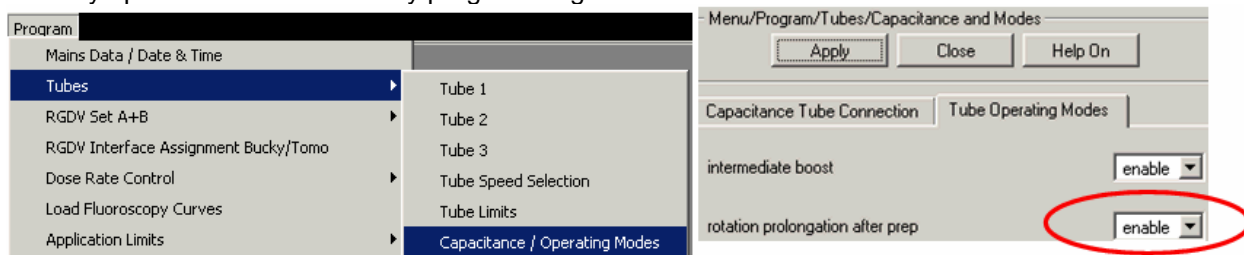

Then the anode keeps rotating for 30 seconds after the last let go of the prep switch and can be repeated as often as necessary.

Recommended especially for children's and A&E rooms.

RF systems require the same setting:

If customers want to make exposure series with individual single exposures they should release fluoro once to keep the anode speed for 30 seconds, otherwise 00XT +10LA might occur again.

Important remark:

The above measures don't work if customers push both steps PREP **and** EXP together.

The exposure request lets the generator behave as if "rotation prolongation after prep" is not enabled.